
WATER PERMEABLE SURFACE

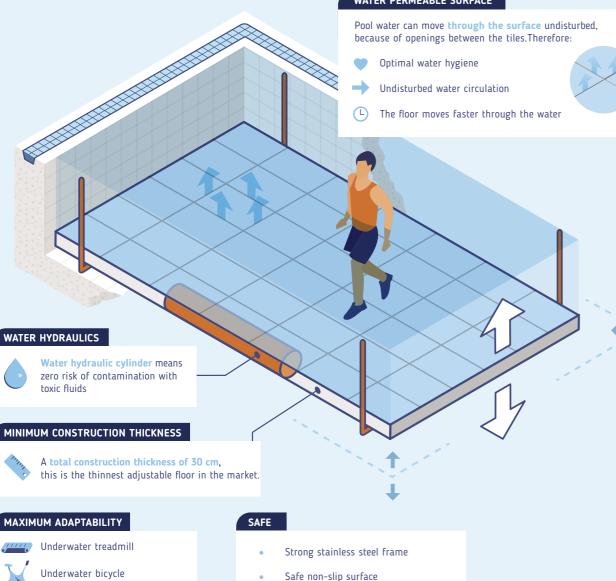

- Water permeable over entire surface
- Safe water hydraulics

DYE TEST

1

Only adjustable floor that passes the EN15288-2:2008 annex A Dye test

Ready for use

Single point

of contact

No toxic

fluids

in 1 week

Additional advantages

Optimal water

maintenance

1000⁴

projects worldwide

hvgiene

Easy

QUALIFICATIONS

11

 \bigcirc

Certified as Class I Medical Device

Obstacle course

Counter current

Underwater massage

Self adjusting steps

Camera's

ISO 9001 Certified manufacturer

ISO 13485 Qualified supplier

RETRO FITTING POSSIBLE

Can be installed in new and existing pools

EWAC Medical Galileistraat 20 1704 SE Heerhugowaard The Netherlands



+31 226 313 457 medical@ewac.nl www.ewacmedical.com

25+

years

lifespan

80% recyclabl

Subject to alterations and errors. No rights may be derived from this brochure

MOVABLE POOL FLOOR



9

www.ewacmedical.com



MOVABLE SWIMMING POOL FLOOR



We believe in aquatic rehabilitation, a therapist should create the optimal exercise environment for the patient, water depth being one of the most important variables.

This is why we have created a movable swimming pool floor that is specially developed for medical applications. Our design having unprecedented water permeability, our adjustable pool floor offers the very best hygienic performance in the market. Very important with vulnerable patients.

In aquatic therapy, you should choose the depth to suit the patient. Age, height and progress among other factors, play an important role in determining the right exercise depth. With our adjustable pool floor, the depth can be changed at the push of a button. The floor can be raised quickly to deck level, so that patients can easily be transferred in- and out of the pool.

Versatile

The EWAC Medical movable pool floor is the most versatile product on the market. We can fit it to any pool, any size. Even in existing pools a adjustable floor can be retrofitted. It can be equipped with a Pooltrack® integrated underwater treadmill, facilitating stationary walking exercises. We offer a wide range of pool accessories for walking and cycling

In new pools, camera's and underwater lighting can be fitted as well as jet streams, counter currents and underwater massage systems, greatly enlarging the treatment possibilities for the therapist.

exercises and water based balance exercises

With our versatile Modular Pool, we can even install the pool quickly without any heavy construction work on your site! This means you will have a single point of contact for the entire hydrotherapy solution.

Hygienic


Our adjustable floor system meets the strictest regulations for water hygiene. In fact it is the only system on the market that meets the requirements of the European EN 15288:2008, annex A dye test. In effect this means that the presence of the movable floor does not affect the water circulation in the pool, as is necessary for proper water disinfection.

Safe

Our adjustable floor has a strong stainless steel AISI 316 L frame, with a safe non-slip polyester surface. Water can move through the entire surface at all times, evenly spreading the circulated water over the pool.

Safe water hydraulics

Movement of the floor is achieved with safe water hydraulics. Compared with oil hydraulic systems, this offers unprecedented protection against contamination with toxic hydraulic fluids in the pool water.

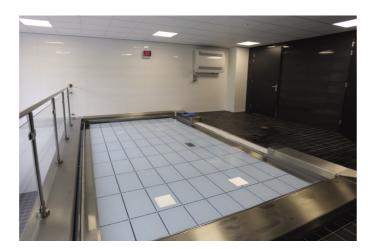

Easy maintenance

All the mechanical parts of the system are easily accessible from the surface of the floor. Service and cleaning can be done by just removing a few of the panels. Contrary to other products in the market, there is no need to hire a technical diving team to perform maintenance underneath the adjustable floor.



EAS

Environmentally friendly

We use more than 80% recyclable materials. At the same time, a adjustable floor helps preserve energy: if it is set to deck level at night, it prevents water evaporation and heat radiation. Therefore it prevents heat loss and reduces CO2 emissions.



Accessories:

The adjustable floor has a range of useful accessories that can accelerate the patient's rehabilitation and give the therapist options to enhance the rehabilitation progress and fun.




Exercise accessories:

- Pooltrack[®] Underwater treadmill
- Underwater bicycle
- Obstacle course
- Removable parallel bars
- Hand rails



Training accessories:

- Aqua fins
- Aqua dumbbells
- Agua noodles

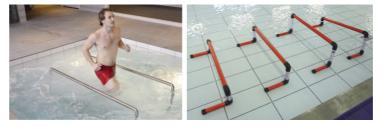







Access options:

- Self-adapting stairs (steps adjusting to the depth of the adjustable floor)
- Gradual entry ramp (for easy entry and exitf o r wheelchairs)



- Patient handling accessories:
- Ceiling mounted hoist
- Mobile chair
- Stretcher / Trolley
- Pool lift





EWAC Medical is Lloyd's Register Quality certified according to ISO 9001 and ISO 13458. EWAC Medical has the privilege to design, manufacture and test Class I Medical products according to the European Medical Device Directive.

- EWAC Medical is certified according to: ISO 9001 ISO 13485
- by Lloyd's Register Nederland



www.ewacmedical.com