



TRANSFER EQUIPMENT

ALWAYS THE BEST CHOICE...

SAFE AND COMFORTABLE

EWAC

**MEDICAL**

*We get you moving*



[www.ewacmedical.com](http://www.ewacmedical.com)



## Universal Ceiling Hoist (UCH)

The universal ceiling hoist is a rail mounted lifting module that is user-friendly and powerful enough to perform lifting and moving procedures in almost all professional care environments. The lifting module is easy to operate and helps to make lifting and moving procedures safe, simple and comfortable for both the user and the nursing staff.

The lifting module charges itself via the rails, regardless of where in the system the lifting module is located. This ensures that the ceiling hoist is always ready for use. This lifting module is specially designed for the working conditions typical of professional healthcare environments and is dimensioned to handle heavy loads on a regular basis. The UCH is the ideal choice for many different lifting and moving procedures, and can also be used to assist in other care tasks, mobilization exercises, training and rehabilitation programs.



Several versions of the UCH lifting module are available, and the module can be configured and adapted to individual needs with regards to functionality, specifications and budget.

The UCH is available in models with a lifting capacity of 200 and 250 kg, respectively. The UCH can be equipped with a built-in motor for horizontal movement along the traverse rail, with the lifting module being operated by hand control or infrared (IR) remote control.



## Universal Sling

The Basic sling is a general lifting sling that is suitable for a wide variety of users and lifting activities. The sling supports the whole body except the head.

The sling is ideal for lifting, transferring and positioning users who can control their head but have less control over their upper or lower body. It can be used from a sitting or lying position - for example, from a bed or wheelchair - The Basic Basic sling is designed as a general lifting sling for most lifting situations.

The Basic Basic sling is available in different sizes and two different materials, polyester or polyester net, which can be used as a bath sling. Both materials are very durable and easy to work with. They allow air to circulate around the body and provide excellent comfort and support for the user.




The standard lifting capacity is up to 255 kg; however, the lifting capacity of the 3XL is 500 kg. On the sling, you will find a label with a QR code that takes you to a mobile site with links to instructional video and user manual.

The sling can be used with all ceiling hoists and mobile elevators.

For more information on how to choose the right sling and size, please refer to the size chart below.

Using the correct sling is essential for efficient transfers to ensure comfort and safety for the user and caregiver.

### Guidelines on choosing the size of a sling

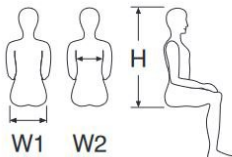ABC Slings	Kids	Kids	Kids	XS	S	M	L	XL	2XL	3XL
	4 - 6	6 - 10	10 - 14							
	H	45 - 60	55 - 70	68 - 80	75 - 82	80 - 87	85 - 92	90 - 97	95 - 102	95 - 102
	W1	25 - 30	28 - 35	33 - 38	33 - 37	36 - 40	39 - 43	42 - 45	45 - 49	48 - 52
	W2 <sup>x)</sup>		60 - 95	65 - 100	70 - 105	75 - 110	80 - 120	90 - 130	100 - 140	110 - 150
	W2 <sup>xx)</sup>		28 - 35	33 - 38						

**W2 = Chest measurements**

x) Guideline for choosing Active Trainer

xx) Guideline for choosing Active Vest Kids

*All dimensions in cm*



## Swimming pool hoisting chair self-propelled


- Seat and back HT padding.
- Lifting eyes on back rest and front legs.
- Separately foldable armrests, ergonomically shaped.
- Armrests in aluminium, white coated.
- Stainless steel swivel castors Ø 125 mm with double brake and plastic 24" castors with brake and hoop.
- Adjustable and foldable footrests.
- Only available in seat height 490 mm.
- Only available in stainless steel 316 polished version available.
- Customization not possible.



## Swimming pool hoisting chair

- Seat and back HT upholstery.
- Lifting eyes on back rest and front legs.
- Separate folding armrests, ergonomically shaped, aluminium, with a white coating.
- Four stainless steel swivel castors, Ø 125 mm with double brake.
- Adjustable and foldable footrests.
- Only available in seat height 490 mm.
- Only available in stainless steel 316 polished
- Customization not possible.



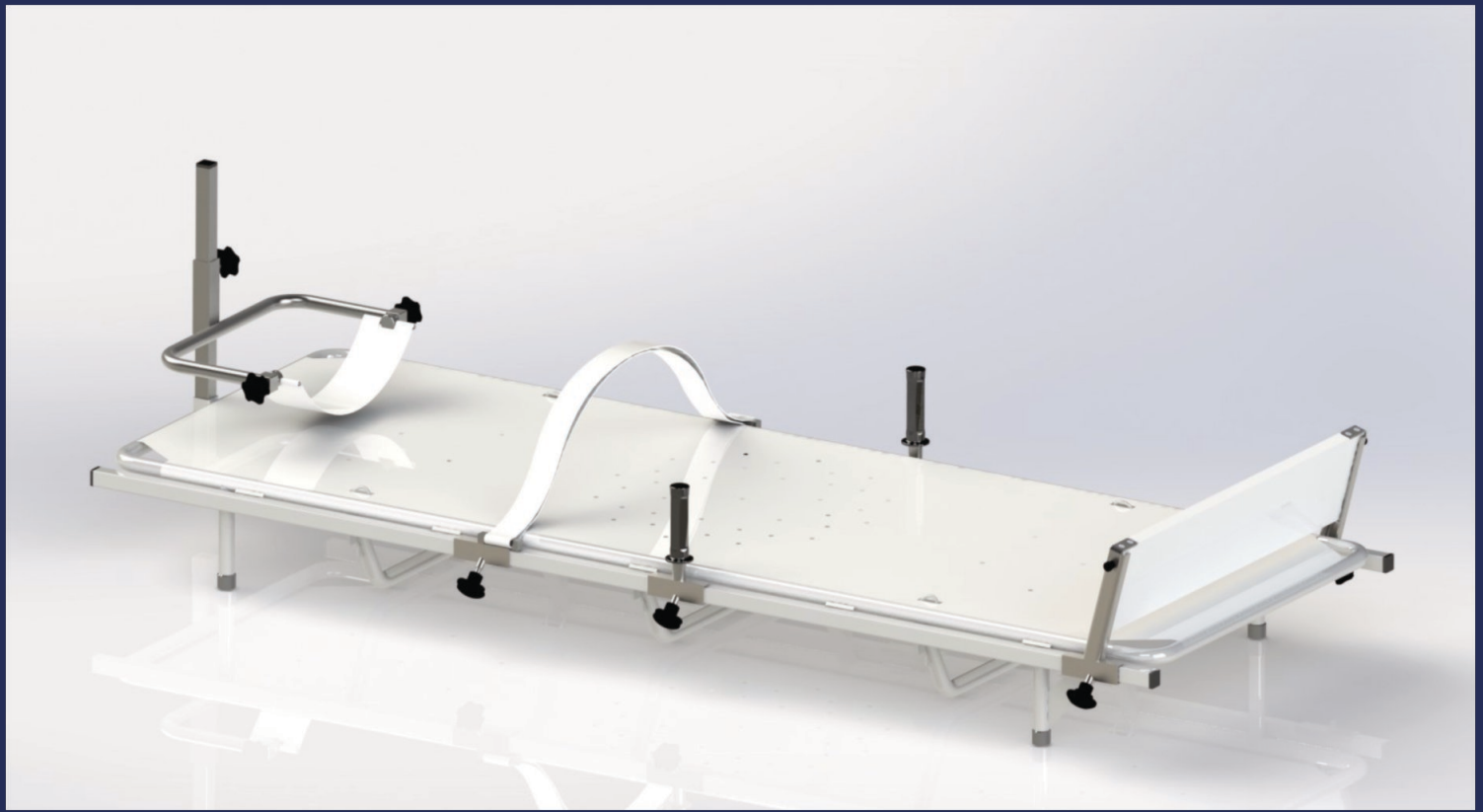

### Stretcher Universal with Trolley

- Lying surface HT padding.
- Four stainless steel swivel castors, Ø 125 mm with double brake.
- Adjustable head and foot rest.
- Set of handles and fixation strap.
- Upper frame is detachable from lower frame and equipped with lifting eyes.
- Only available in the indicated sizes with length 1905 mm.
- Only available in stainless steel 316 polished version.
- Maximum user weight 150 kg.
- Weight: 11 kg.
- Dimensions: 1905 x 660 x 250 mm.



## (Combi) bath lifter

Multifunctional bath seat for optimal support.

In a care bathroom a bath lifter is indispensable! A bath lifter ensures that people with reduced mobility can continue to use the bath. With the bath lifter, clients can get in and out of a high/low bath without any effort. This is important in a (nursing) care facility, because high/low baths are not only used to care for clients, but they also offer a bit of relaxation and therapy.

### Ergonomic bath lift

Our bath lifter is excellent for safe and comfortable transfer from bed to bath. Because the bath lifter is electrically adjustable in height, the caregiver can set the right working height. The bath lifter is available with a bath seat and/or bath stretcher.

During bathing, the client does not need to leave the chair or stretcher. This increases the efficiency for both client and caregiver.



### Features bath lifter:

- Seat and lying surface of the bath elevator are made of soft, durable and waterproof material.
- Seat has a toilet opening with a lid.
- With two removable sides, the seat of the bath elevator becomes a stretcher where the client can lie on.
- New compact chassis with good mobility.
- Innovative swivel castors.
- Comfortable backrest.





## Pool Hoist

The pool hoist is specially designed for use at the edge of a pool. This hoist system makes it possible for people to enter the water who would not normally be able to get into the pool via stairs or a ladder. The pool hoist helps with this in an easy and safe way.



## SPECIFICATIONS

- For persons weighing up to 150 kg.
- Shock-free height adjustment by electric drive.
- Integrated electronically controlled battery charger.
- Mobile pneumatic hand control with a 2 m long control cable so that the therapist can always be close to the client - even in the water.

## USES

- The Pool Hoist can be equipped with a seat or bench and with or without a trolley.
- Mounted with 4 composite anchors on the floor or balustrade.
- Alternatively, the removable Pool Hoist can be mounted with a ground socket.
- Individual adaptation to different types of pools, with or without balustrade.
- Robust due to stainless steel coated with plastic.







EWAC Medical is certified according to:  
ISO 9001  
ISO 13485  
By Lloyd's Register The Netherlands

EWAC Medical  
Galileistraat 20  
1704 SE Heerhugowaard  
The Netherlands



+31 (0)226 313 457  
medical@ewac.com  
www.ewacmedical.com

Subject to alterations and errors. No rights may be derived from this brochure.